


STRUCTURAL FLOOR PLAN

SCALE: 1/4" = 1'-0"

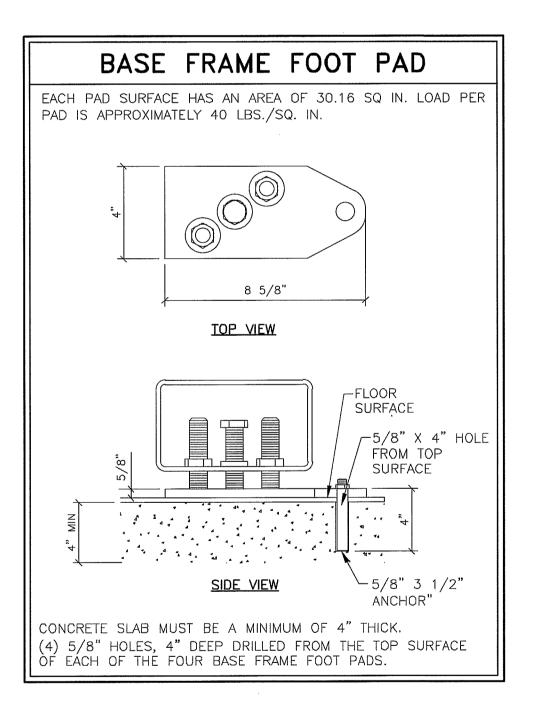

1) THE CUSTOMER/CONTRACTOR SHALL FURNISH AND INSTALL ALL STRUCTURAL SUPPORT MEMBERS AND NEEDED HARDWARE FOR THE INSTALLATION OF THE SIEMENS EQUIPMENT.

STRUCTURAL NOTES

2) THE OVERHEAD STRUCTURAL SUPPORT SYSTEM SHALL BE FIXED, RÍGID AND BRACED FOR SWAY.

3) ALL STRUCTURAL SUPPORT MEMBERS SHALL BE TRUE, SQUARE, LÉVEL. PARALLEL AND COPLANAR WITH RESPECT TO EACH OTHER, WITH A HORIZONTAL STRUCTURAL SUPPORT MEMBER TO BE LOCATED AND SET WITH A TRANSIT.

4) ALL STRUCTURAL SUPPORT DETAILS SHOWN ARE SAMPLE DETAILS BÁSED UPON TYPICAL AND STANDARD BUILDING PRACTICES AND ARE NOT INTENDED AS ACTUAL CONSTRUCTION DETAILS. ALL CONSTRUCTION DETAILS AND SUPPORT CALCULATIONS SHALL BE PREPARED BY A PROFESSIONAL STRUCTURAL ENGINEER AT THE CUSTOMER'S EXPENSE. IN THE EVENT AN EXISTING SUPPORT SYSTEM IS TO BE USED, IT WILL BE THE CUSTOMER'S RESPONSIBILITY TO VERIFY THE INTEGRITY OF


5) WHERE SHOWN ON THE 1/4" STRUCTURAL FLOOR PLAN, THERE ARE ON OCCASION MOUNTING FRAMES FURNISHED BY SMS. THESE FRAMES ARE TO BE SET BY THE CONTRACTOR, UNDER THE SUPERVISION OF SMS PERSONNEL. THE CUSTOMER/CONTRACTOR IS RESPONSIBLE FOR ALL FRAMES INSTALLED BY HIM TO BE WATER LEVEL AND ANCHORED


6) ALL CEILING FIXTURES (I.E. AIR SUPPLY GRILLES, AIR RETURN GRILLES, EXHAUST GRILLES, SPRINKLER HEADS, INCANDESCENT AND FLUORESCENT LIGHT FIXTURES, INTERCOM SPEAKERS, MEDICAL GAS COLUMNS, ETC.) SHALL BE INSTALLED FLUSH MOUNTED WITH THE FINISHED CEILING TO PROVIDE FREE AND UNRESTRICTED TRAVEL OF THE SMS CEILING MOUNTED EQUIPMENT.


7) THE BOTTOM SIDE OF THE UNISTRUT CEILING GRID AND ANY CEILING MOUNTED SUPPORT PLATES ARE TO BE INSTALLED FLUSH WITH THE FINISHED CEILING. THE CUSTOMER/CONTRACTOR SHALL ALSO PROVIDE COVERSTRIPS FOR THE UNISTRUT.


8) THE STRUCTURAL PLANNING AS SHOWN ON THE 1/4" STRUCTURAL PLAN HAS BEEN COORDINATED WITH THE EQUIPMENT LOCATION AS SHOWN ON THE 1/4" EQUIPMENT LAYOUT PLAN. FOR THIS REASON, ANY DEVIATIONS FROM THE STRUCTURAL PLANNING AS SHOWN MUST BE APPROVED BY SMS PLANNING DEPARTMENT.


9) THE STRUCTURAL ENGINEER OF RECORD SHALL BE RESPONSIBLE FOR THE DESIGN AND DETAIL OF FLOOR, WALL AND CEILING STRUCTURES IN ACCORDANCE WITH THE WEIGHTS, MOMENTS AND FORCES AS SHOWN ON OUR STRUCTURAL CALCULATIONS, OR INFORMATION, IN CONSIDERATION OF FORCES AS DETERMINED PER LOCAL GOVERNING BUILDING CODES.

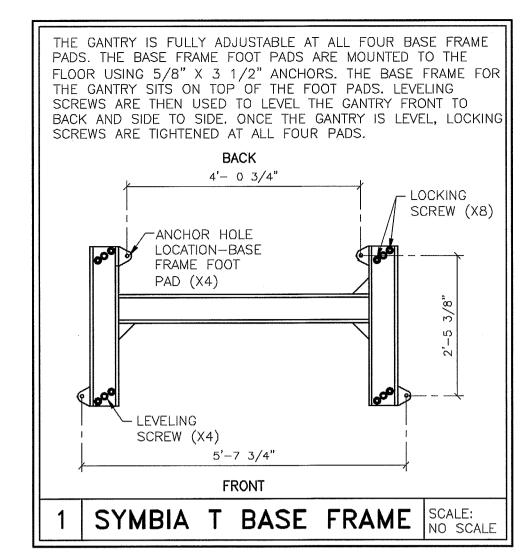

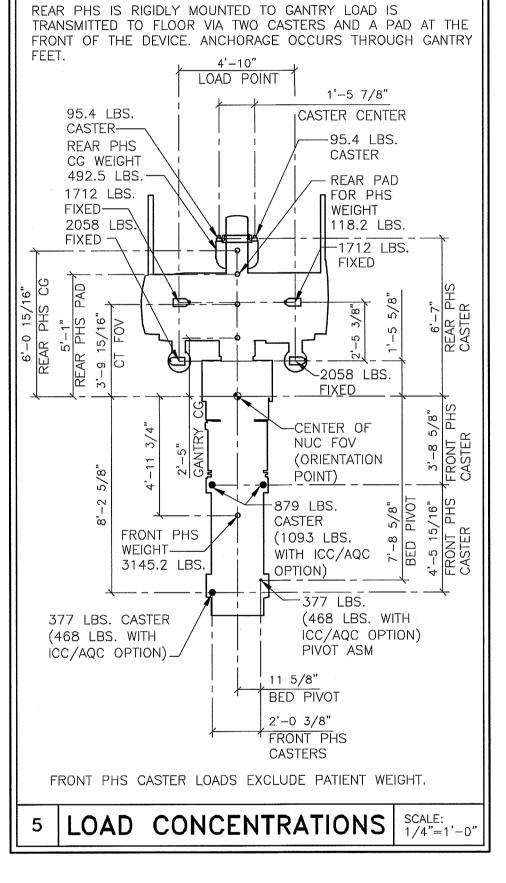

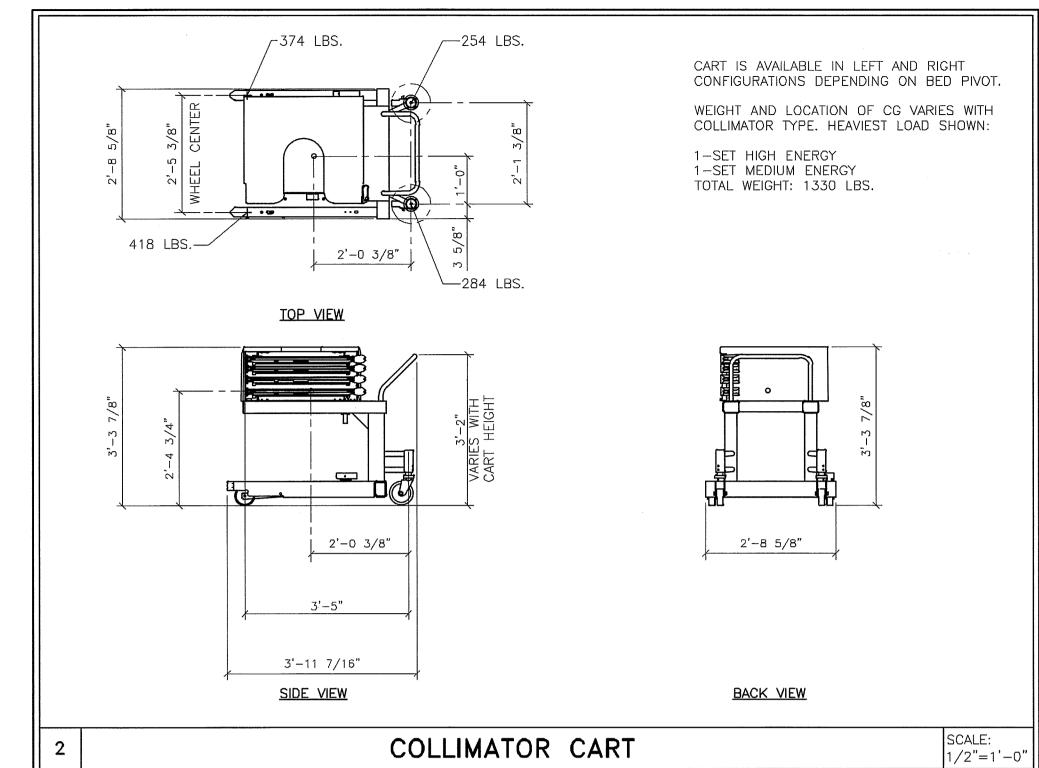

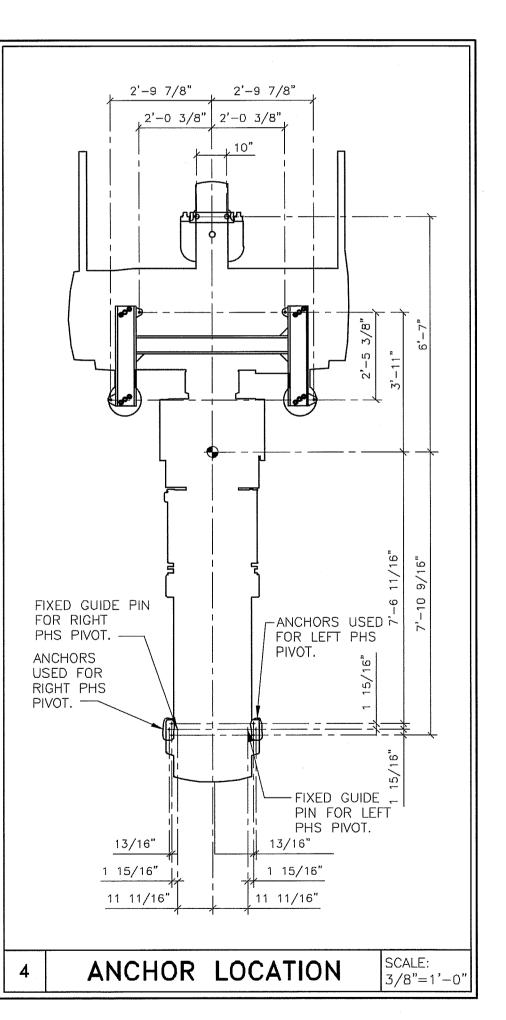

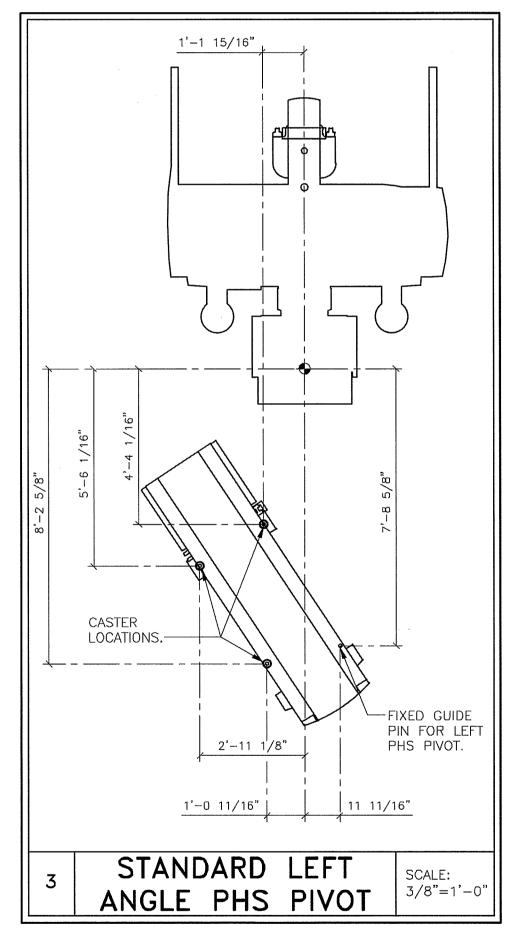

REV 2

MUST BE PLACED AS TO NOT COLLIDE WITH WARNING LIGHT.

FINISHED ROOM HEIGHT SYMBIA T, T2, T6 OR T16 MINIMUM 8'-0" MINIMUM 8'-0" SYMBIA T, T2, T6 OR T16 WITH CEILING MOUNTED COMPONENT MAXIMUM 12'-0" OTHER THAN RADIATION ON LAMP CONSIDER THE WARNING LIGHT WILL BE PLACED ON TOP OF THE PATIENT BOOM. ANY OTHER CEILING MOUNTED COMPONENT

PROJECT MANAGER: RICH DEISTER
TEL: (207) 712-3205 **SIEMENS** VMAIL: (207) 929-3776 EMAIL: rich.deister@siemens.com MAINE MEDICAL CENTER 22 BRAMHALL STREET, PORTLAND, ME 04102 NEW SYMBIA SUITE - SYMBIA T 02/10/14 REVISED ELECTRICAL AND REMOVED THE USE OR REPRODUCTION OF PROJECT #: SYMBIA T PROCESSING WORKPLACE. THIS TITLE BLOCK WITHOUT SIEMENS AUTHORIZATION WILL 1303838 R-101R(A) VERSION DATED 12/ RESULT IN PROSECUTION UNDER APPROVED BY THE CUSTOMER FOR FIN FULL EXTENT OF THE LAW. ALL RIGHTS ARE RESERVED. DESCRIPTION M. COPLEY AS NOTED REF. #: **400**-347620 -ISSUE BLOCK-01/14/14

ATTENTION:

- THIS DRAWING IS DESIGNED TO CONFORM TO FEATURES AND EQUIPMENT REQUIREMENTS PRESENTED AT THE TIME OF THEIR PREPARATION. SINCE BOTH THESE FACTORS ARE SUBJECT TO DESIGN MODIFICATION, THEY ARE NOT TO BE USED FOR CONSTRUCTION PURPOSES. - THIS SET OF PLANS REPRESENTS A COMPLETE SET OF DETAILS AND SHOULD NOT BE SEPARATED.

-IT IS RECOMMENDED THAT THE SIEMENS DRAWINGS BE INCORPORATED WITH THE CONSTRUCTION DOCUMENTS FOR REFERENCE.

- ALL DIMENSIONS SHOWN ON THIS DRAWING ARE FROM FINISHED SURFACES.

- THIS DRAWING DOES NOT PROVIDE RADIATION SHIELDING REQUIREMENTS FOR X-RAY AND ASSOCIATED EQUIPMENT. THE CUSTOMER IS RESPONSIBLE FOR CONSULTING WITH A REGISTERED RADIATION PHYSICIST TO SPECIFY RADIATION PROTECTION.