

PELTON & CRANE RF42 CABINET UTILITIES

— 1/2" FEMALE ADAPTOR

CONTRACTOR.

PROVIDED & INSTALLED BY

- AIR LINE BY CONTRACTOR

ANCHOR TERMINATION

FITTINGS PROVIDED AND

OR A CROSS BRACE.

FITTINGS SECURELY TO A STUD

INSTALLED BY CONTRACTOR.

SEE ACTUAL VALVE TO ADJUST THE MOUNTING

SHEETROCK.

POSITION OF THE 1/2" FPT

ADAPTOR SO THE TRIM RING IS

SECURED FIRMLY AGAINST THE

AIR VALVE PROVIDED AND

3/8" COMPRESSION FITTING

DISCONNECT PROVIDED &

OR OPTIONAL FEMALE QUICK

INSTALLED BY CONTRACTOR. —

INSTALLED BY CONTRACTOR.

 $\overline{}$ HIS DETAIL IS FOR INFORMATIONAL PURPOSES ONLY. DENTAL DEALER WILL PROVIDE INSTALLATION TEMPLATES FOR ANY EQUIPMENT IN QUESTION. REFER TO MANUFACTURER'S INSTALLATION GUIDE/TEMPLATE FOR ADDITIONAL INFORMATION. -LICENSED PLUMBERS AND ELECTRICIANS SHALL INSTALL ALL UTILITIES.

-ALL UTILITIES MUST CONFORM WITH PREVAILING LOCAL CODES. NEITHER DENTAL DEALER NOR MANUFACTURER CAN BE HELD RESPONSIBLE FOR AN INSTALLATION WHICH DOES NOT MEET LOCAL ELECTRICAL AND PLUMBING CODES.

DENTAL STERILIZER (STATIM)

(D-15B) KAVO QUATTROcare

-VERIFY WITH EQUIPMENT REP WHICH UTILITIES WILL BE RUN SUB-FLOOR VS. IN WALL.

AIR: PLUMBING CONTRACTOR TO PROVIDE AND INSTALL A 1/2" COPPER (OR ACCEPTED) SUPPLY LINE FROM AIR COMPRESSOR TO THIS LOCATION. STUB-UP 1-2" FROM FLOOR VACUUM SUPPLIER). REDUCE TO A 1/2" SCH. 40 PVC (OR AND TERMINATE WITH A 90° 1/2" x 3/8" COMPRESSION SHUT-OFF VALVE. FLUSH SYSTEM CLEAN PRIOR TO EQUIPMENT FLOOR FOR CONNECTION TO EQUIPMENT. INSTALLATION.

NOTE: MAXIMUM HEIGHT OF AIR VALVE SHOULD NOT EXCEED 5" (WHEN FULLY OPENED).

ELECTRICAL: ELECTRICAL CONTRACTOR TO PROVIDE AND INSTALL A 1/2" FLEX CONDUIT AND BOX WITH 115 QUAD (OR EQUAL) RECEPTACLE. WIRE BOX PER CODE WITH TOP NO HIGHER THAN 4-1/2" A.F.F.

CENTRAL VACUUM: PLUMBING CONTRACTOR TO PROVIDE A BRANCH LINE TO LOCATION SHOWN (SIZE PER CENTRAL COPPER) RISER AT FLOOR. TERMINATE 1-2" ABOVE FINISHED

CHASE: ELECTRICAL CONTRACTOR TO PROVIDE A GRAY ELECTRICAL PVC PIPE CHASE WITH PULL-STRING SUB-FLOOR (SEE SPECS 4, 4A SHEET SE.1 FOR SIZE AND LOCATIONS). 6" MIN. RADIUS FOR ALL BENDS. STUB-UP FLUSH WITH SUB-FLOOR, MARK LOCATION OF PIPE AT FINISH.

OUTLET LOCATIONS ARE APPROXIMATE

AND SHOULD BE VERIFIED BY HENRY

DENTAL STERILIZER PROVIDED AND

INSTALLED BY HENRY SCHEIN DENTAL

ELECTRICAL: ELECTRICAL CONTRACTOR

TO PROVIDE AND INSTALL A BOX WITH A

DEDICATED 115V, 15A CIRCUIT & SINGLE

PROVIDE A 9/32" HOLE IN COUNTER

SLIDE-OUT WATER RESISTANT SURFACE

MUST BE CAPABLE OF A 100 LB LOAD.

- CONDENSING TANK (PROVIDED WITH

AIR: PLUMBING CONTRACTOR TO PROVIDE AND INSTALL AN AIR SUPPLY LINE AT 44" A.F.F. STUB-OUT FROM WALL AND TERMINATE WITH A 90° 1/4"

ELECTRICAL: ELECTRICAL CONTRACTOR

120V, 20A CIRCUIT & SINGLE RECEPTACLE

TO PROVIDE AND INSTALL A BOX WITH

OUTLET LOCATIONS ARE APPROXIMATE

AND SHOULD BE VERIFIED BY HENRY

QUATTROcare UNIT PROVIDED AND

INSTALLED BY HENRY SCHEIN DENTAL.

MIN. 8" SIDE CLEARANCE FROM ALL

20" SIDE CLEARANCE FROM ALL IGNITION SOURCES, LIGHT SWITCHES,

UNIT 11" DEEP WHEN CLOSED.

WALL AND CABINET SURFACES. MIN.

PROVIDE AN ADDITIONAL 8" IN FRONT OF UNIT FOR WHEN FLAP IS OPEN.

SCHEIN EQUIPMENT REP.

AND MAIN OUTLETS.


COMPRESSION STOP.

AT 44" A.F.F.


RECEPTACLE AT 30" A.F.F.

BEHIND UNIT TO RUN TUBING.


SCHEIN EQUIPMENT REP AND CABINET

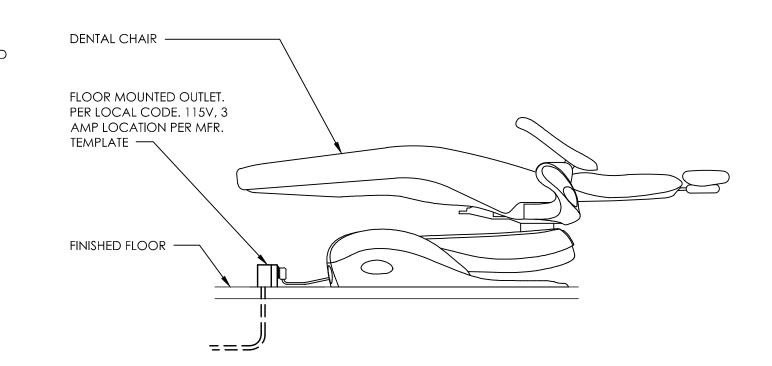


FLOOR-MOUNTED CHAIR OUTLET NOT TO SCALE

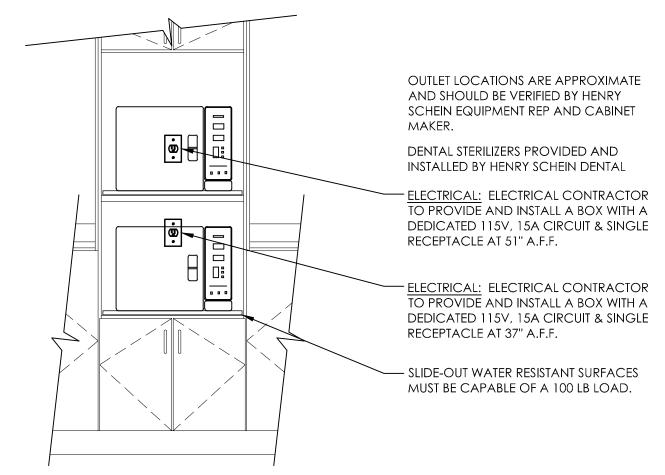

STACKED STERILIZERS (AUTOCLAVES)

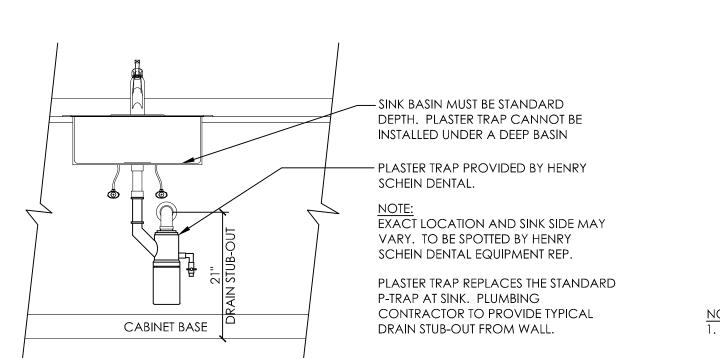




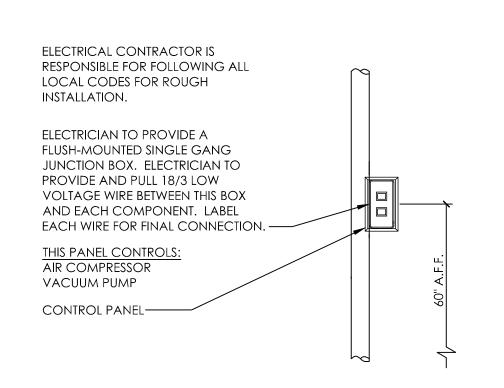


CONTROL PANEL PROVIDED BY HENRY-SCHEIN DENTAL, FINAL CONNECTION BY HENRY SCHEIN DENTAL. SEE MANUFACTURER'S INSTALLATION MANUAL FOR MORE DETAIL.




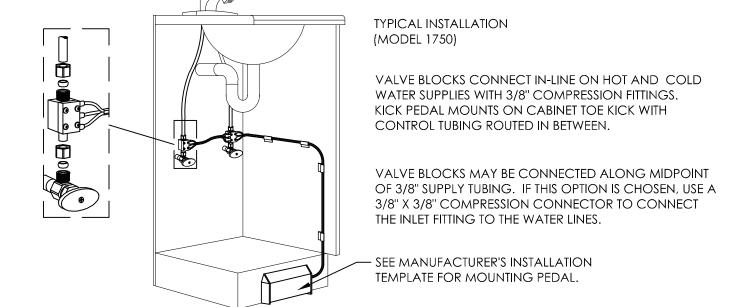




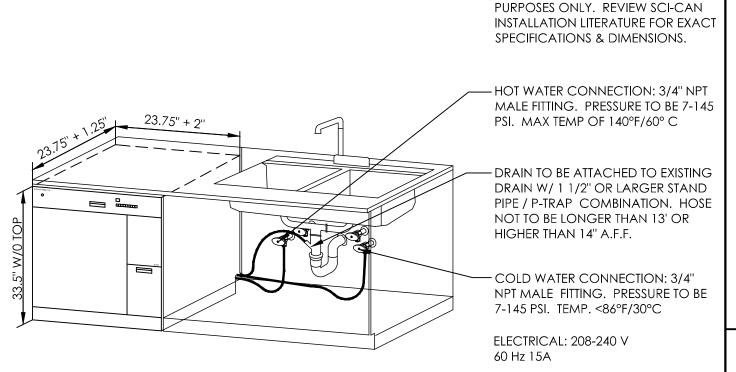

hydrim instrument washer

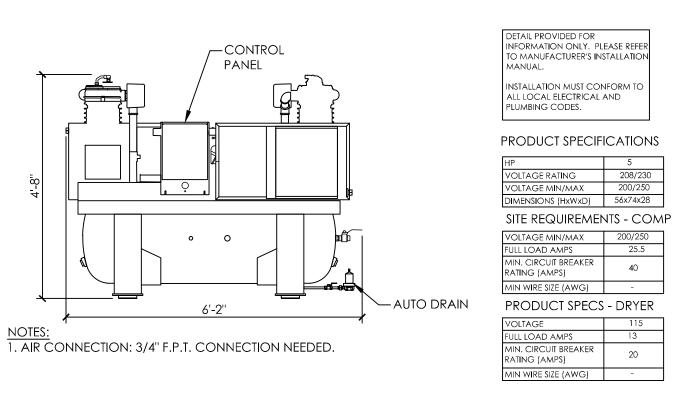



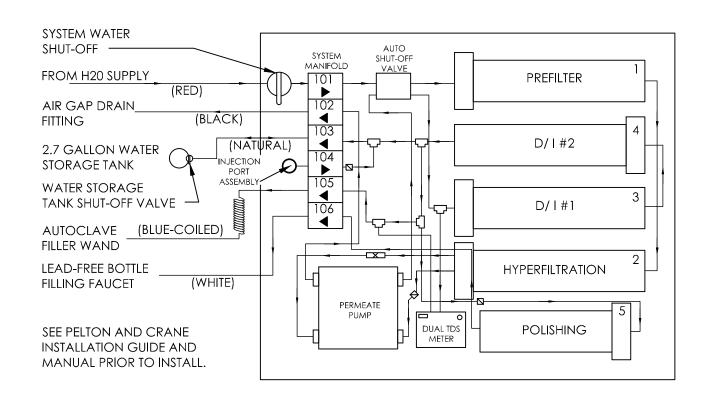

*IF NOT ALREADY MOUNTED IN ANOTHER PRODUCT (STERI CENTER) MOUNT NEAR COLD WATER SOURCE AND 1-1/2" DRAIN.

PELTON & CRANE VISTAPURE

*ALWAYS FOLLOW LOCAL PLUMBING CODES.



THIS DRAWING FOR ILLUSTRATION

CHODON


HENRY SCHEIN REP JESSE LEVEILLE CENTER: BOSTON, MA PHONE #: (207) 233-4978


DRAWING NAME: PROSTHO-3d-F PROJECT START DATE: 06/06/2013 | FINALS START DATE: 08/15/2014 DRAWN BY: AKL/CAS FINALS BY:

CHECKED BY: **REVISIONS:** 08/22/14 AKL 10/22/14 AKL INT.SQ.FT.=

SCALE: SHT. SIZE: D SEE DTL. MECHANICAL DETAILS SPE

© 2014 by HENRY SCHEIN, INC.

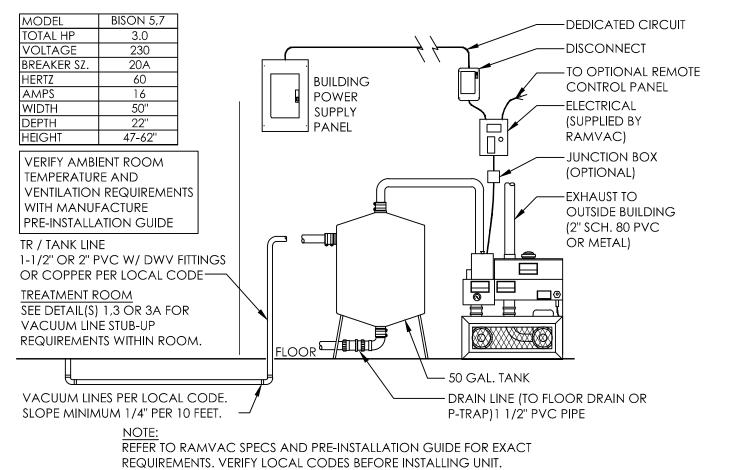

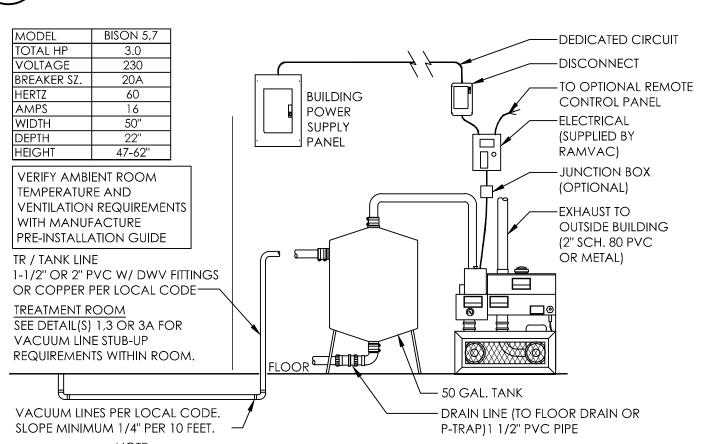

RAMVAC BISON 5 OR 7 VACUUM